

Extractors



Liquid Solid Extraction involves the elution of a solid sample with pure solvent vapour. The sample is held in the Extractor body either directly or in a thimble. The siphoning type of Soxhlet Extractor finds the widest application in laboratory work. Borosil Soxhlet Extraction Apparatus is available in six sizes. The component parts of the apparatus are provided with interchangeable joints for greater flexibility.



3740 - EXTRACTORS Soxhlet, Interchangeable Joint

- Complies with IS 10640
- Made of borosilicate glass
- Designed with a bulb in the siphon tube near the top, to aid the siphoning action
- The vapour tube also protects the siphon tube from accidental damage

Product Code	Capacity ml	Interchangeable Joints		Quantity Per Case	Price / Piece ₹
		Socket	Cone		
3740013	60	34 / 35	24 / 29	5	778.00
3740016	100	40 / 38	24 / 29	5	1100.00
3740020	200	50 / 42	24 / 29	5	1328.00

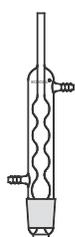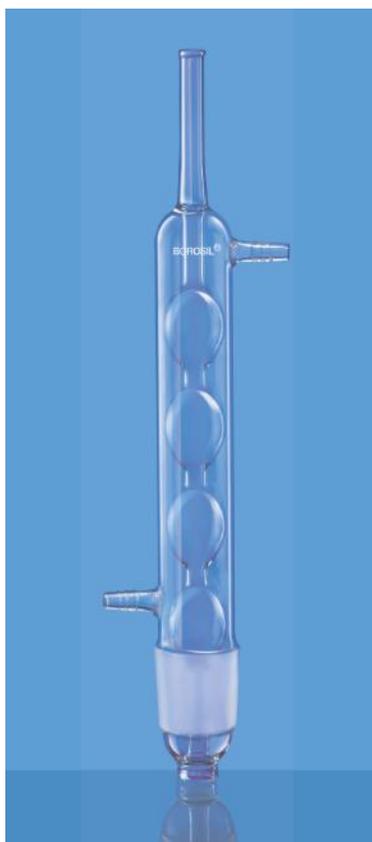

3741 - CONDENSERS, ALLIHN For Soxhlet Extraction Apparatus


- Made of borosilicate glass
- Greater surface area than the Liebig Condenser owing to a series of bulbs, resulting in higher cooling capacity
- Ideal for laboratory scale refluxing

Product Code	Capacity ml	Interchangeable Joints	Quantity Per Case	Price / Piece ₹
		Cone		
3741013	Small	34 / 35	5	637.00
3741016	Medium	40 / 38	5	794.00
3741020	Large	50 / 42	5	1145.00

3840 - EXTRACTION APPARATUS

Soxhlet Complete With Allihn Condenser Interchangeable Joint

- Complies with IS 10640
- Used for separation of the extract from the liquid or solid
- Greater surface area than the Liebig Condenser owing to a series of bulbs, resulting in higher cooling capacity
- This method enables effective separation, with relatively small quantities of solvent
- Ideal for laboratory scale refluxing

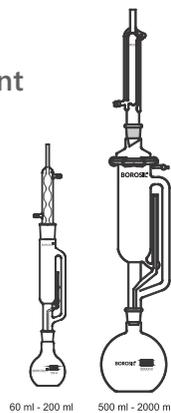

60 ml - 200 ml 500 ml - 2000 ml

Product Code	Capacity ml	Flask Size ml	Quantity Per Case Set	Price / Piece ₹
3840013	60	150	1	1615.00
3840016	100	250	1	1901.00
3840020	200	500	1	3132.00
3840024	500*	2000	1	11362.00
3840029	1000*	3000	1	16972.00
3840030	2000*	5000	1	28652.00

Spares are available on request

For flask refer Cat No. 4100

* Extractor Supplied are with double surface condenser for higher cooling capacity & flange type soxhlet extractors

